FREEDOM
BY ACCESSIBILITY PROFESSIONALS

INSTALLATION INSTRUCTIONS

APHHGBLS Handheld Shower Kit with Glide Bar Installation





Customer Support: 1.877.947.7769
Email: apsupport@accessibilitypro.com

APFHHGBLS Handheld Shower Kit Installation Instructions

Components and General Assembly

1. P/N SEA-SB30-TW-125-SA - 30" Grab Bar; 1.25" outer diameter; stainless steel with satin finish
2. P/N Slider02-CH - Chrome Slide Bracket
3. P/N HS2462CP - Hand Shower; on/off button; chrome finish
4. EPDM Gasket
5. P/N BBS1250x28 - Supply Elbow
6. P/N 4993 - Vacuum Breaker; Brass
7. P/N F14SS60 - Shower Hose, 60" Stainless Steel
8. On/Off Button (flow rate controlled per *on/off* button on body)
9. Mounting Screws
10. P/N NEOPERL 1.8 GPM Flow Restrictor

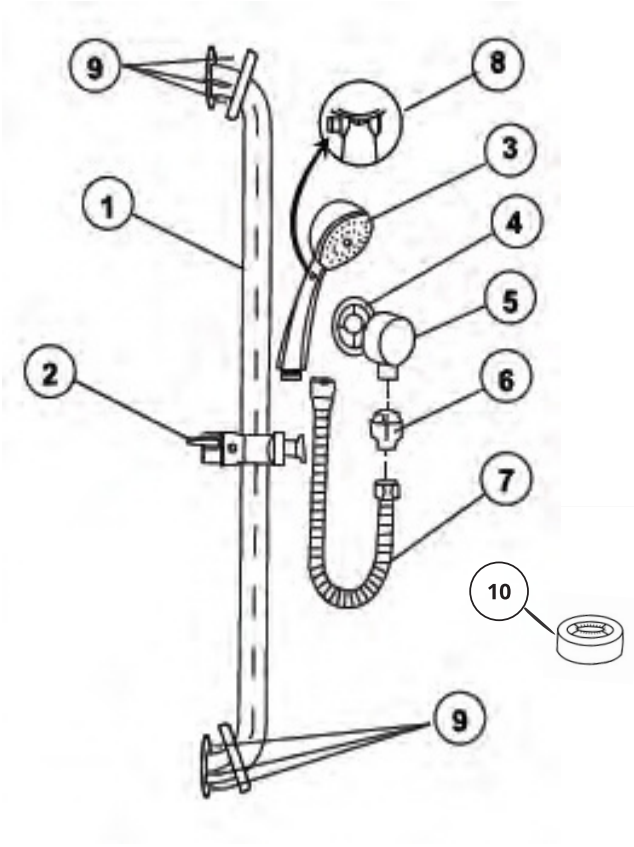

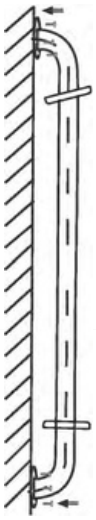











IMPORTANT Have the Grab Bar connected to the wall by a licensed professional.

APFHHGBLS Handheld Shower Kit Installation Instructions

Step 1) Connect the Slide Bar Bracket to the Grab Bar.

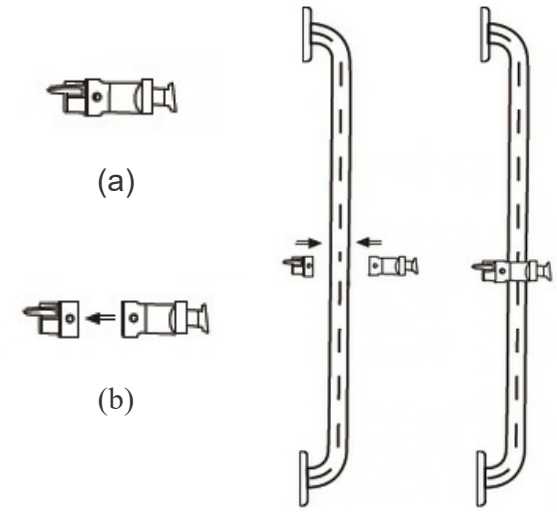

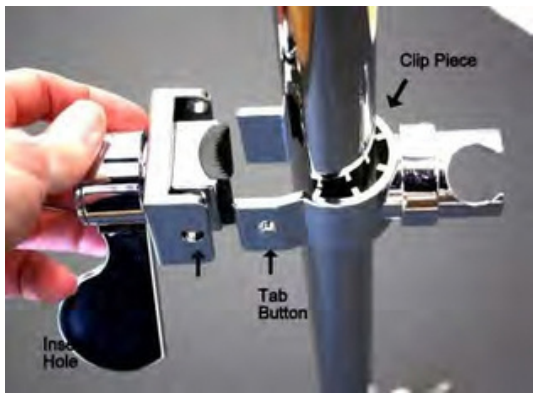


Step 2) The connection to the Water Output Pipe should protrude 10mm-15mm from the shower wall. This connection needs to be wrapped in waterproof tape to ensure a tight seal. Screw the Supply Elbow to the Water Output Pipe.



Step 3) Connect the Supply Elbow to the Vacuum Breaker.






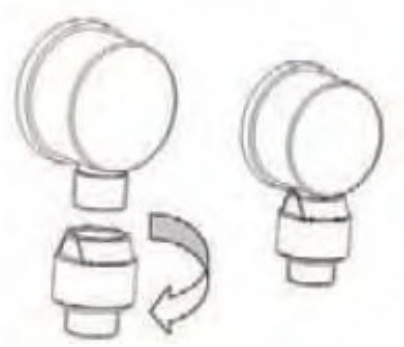

Step 4) Connect the Hose to the Vacuum Breaker. Make sure the rubber washer is placed inside the hose connection.





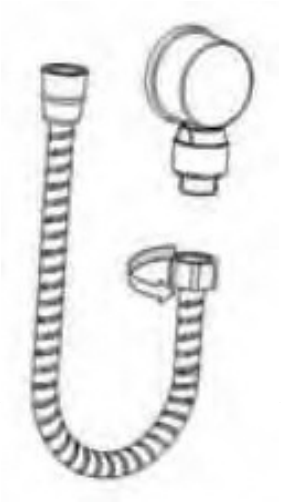


Step 5) Connect the Hand Shower with the Hose as shown. Make sure the rubber washer is placed inside the hose connection if the project does not require flow restriction to 1.8 GPM.

NOTE: If the installation is in California or another location requiring 1.8 GPM max flow, insert the NEOPERL 1.8 flow restrictor into the end of the handshower in place of the rubber washer. Flow restrictor should be inserted blades first so that the green plastic screen is facing outward, as shown below.






Pictured with flow restrictor

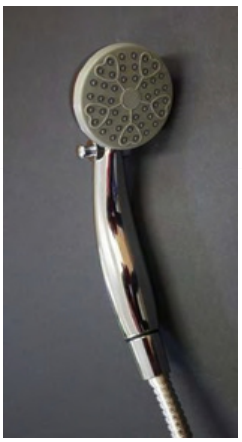


Installing the Sliding Shower Holder


Lever: Moves to loosen and tighten the connection with a Grab Bar in order to slide the attachment up & down
Clip for Attaching to Bar: Slides over a Grab Bar
Clip for holding the handheld shower: Slides over the Handheld Shower
Tab Buttons & Insertion Hole: The tabs click into the insertion holes, locking the two pieces together

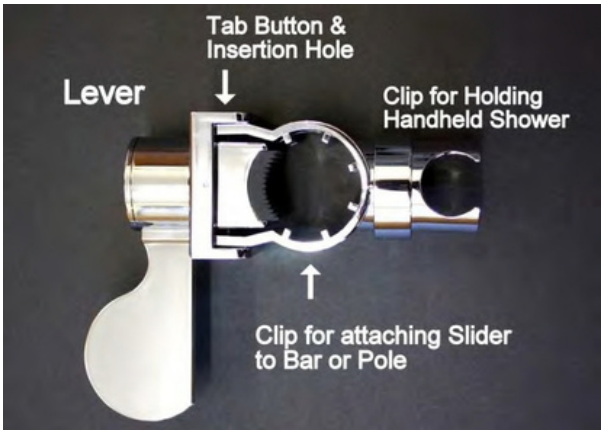


Step 1
Separate the clip piece from the lever piece by depressing both sides of the clip piece at the same time, so the tab buttons will come free of the insertion holes.

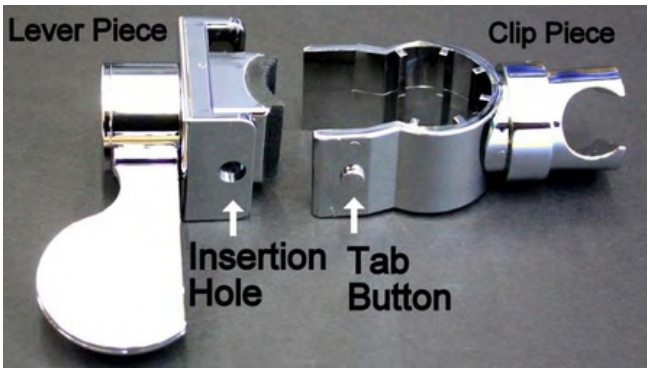


Step 2
Push the arms of the clip piece around the grab bar. The clip piece should fit snuggly around the Grab Bar.

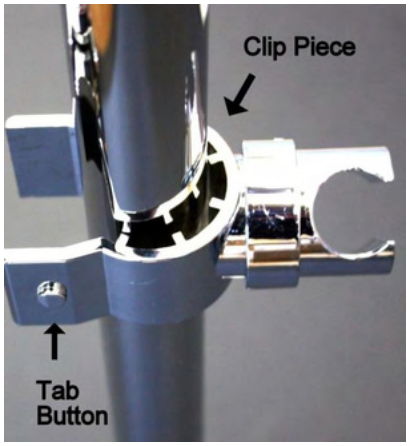

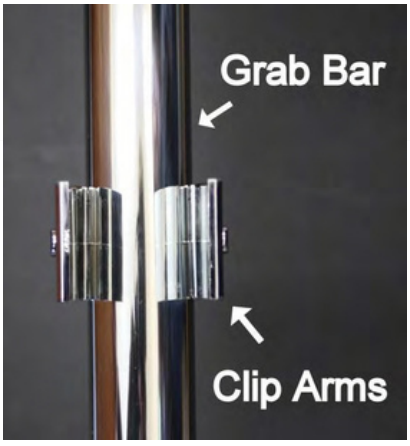


Step 3

Attach the lever piece by sliding the tabs lock into the holes. This will lock the two pieces into place around the Grab Bar.

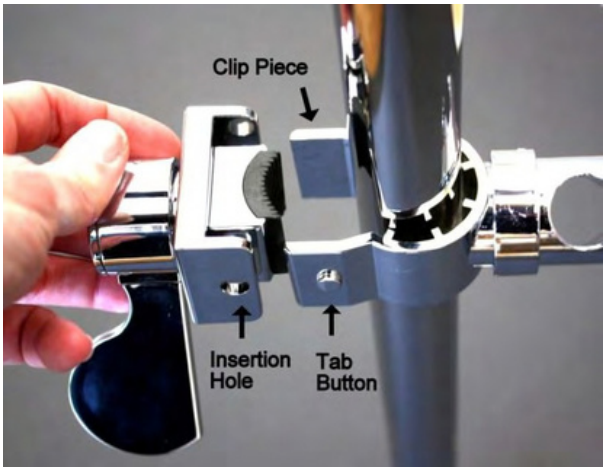

Step 4

When the lever is in the up position, the two piece assembly with slide up or down the Grab Bar, allowing the shower head to be located at the most comfortable height.
Push the lever down to lock the shower into place. Enjoy!

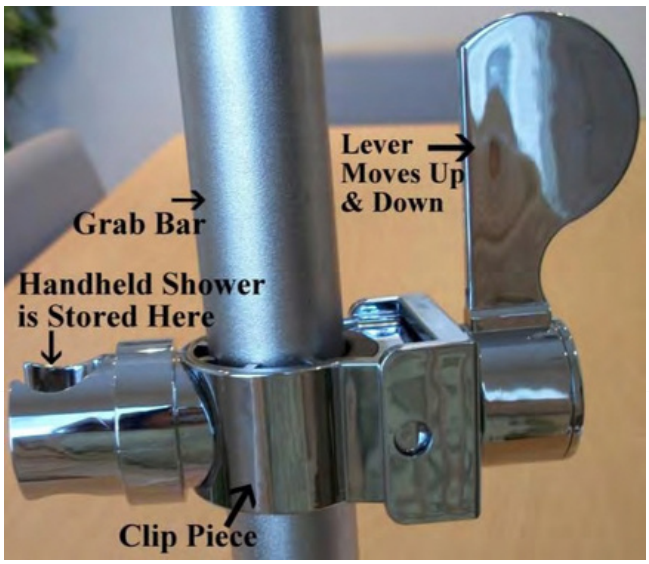

This concludes the installation of the handheld shower kit with grab bar.